



# Microfluidics™

Superior Knowledge | Superior Results



**M815**  
Microfluidizer®  
Pilot-Scale  
Processor Series

# M815 Microfluidizer® Pilot-Scale Processor Series


## For Processing Pilot & Small Production Batches

**Microfluidizer** technology efficiently converts fluid pressure into shear forces, leading industry performance standards in high-pressure processing.

A unique solution to maintaining consistent process pressure ensures 100% of your material gets exactly the same treatment. Whether you are working with small-scale lab batches or production volumes, the Microfluidizer processor is unmatched in submicron particle/droplet size reduction, cell disruption, product yield and guaranteed process scale-up.



High-efficiency cell disruption with minimal protein denaturation

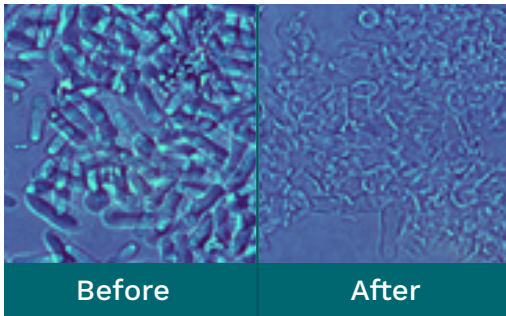

### Recommended for:

- ◆ Emulsions
- ◆ Dispersions
- ◆ Liposomes
- ◆ Cell Disruption
- ◆ Fine Particle Deagglomeration

### Unique Benefits of the M815

- ◆ Produces product flow rates up to 1200 ml/min at 2068 bar (30,000 psi)
- ◆ Has small batch capability; handles a minimum sample size of 1.5 l
- ◆ Features a low product hold-up volume (1 l)
- ◆ Standard with 7" touchscreen HMI that displays process pressure & temperatures
- ◆ Is CIP process capable
- ◆ Integral feed pump
- ◆ Integral heat exchanger
- ◆ Lockable casters, standard door width for easy mobility
- ◆ Meets CE compliance standards
- ◆ Facilitates non-destructive processing of heat-sensitive materials
- ◆ Has cost-effective production capability
- ◆ Assures batch-to-batch process repeatability



# M815 Microfluidizer® Pilot-Scale Processor Series

## Operating Principle

Like all Microfluidizer processors, the M815 utilizes the fixed-geometry Interaction Chamber™ and constant pressure pumping system. This technology allows users to achieve smaller particle sizes — with more uniform distribution and scale-up guaranteed (both from lab-scale and to larger scale units) — than can be obtained with other methods.

The M815 models were designed to bridge the gap between the lab-scale (M110P) and production-scale (M700 Series) models.

## Standard Features

- ◆ Product feed pump, with pressure gauge & purge valve
- ◆ Diamond Interaction Chamber™ (IXC™)
- ◆ Ceramic Auxiliary Processing Module™ (APM™)
- ◆ Ceramic plunger with seal quench for extended seal life
- ◆ Maximum feed temperature of 158°F (70°C)
- ◆ Heat exchanger
- ◆ Stainless steel enclosure
- ◆ Gauges for measuring hydraulic drive pressure, hydraulic oil level & temperature
- ◆ Mounted on locking casters for portability, standard door width

## Specifications

Pressure Range	Up to 30,000 psi (2068 bar)
Flowrate Range	1.0-1.2 l/min
Dimensions (L x W x H)	59" x 34" x 79" (150 cm x 86 cm x 201 cm)
Weight	1950 lbs (886 kg) with oil 1800 lbs (818 kg) without oil

**NOTE:** Some options may change overall dimensions & weight of the machine.

## Options

- ◆ Motor starter panel
- ◆ RTD temperature sensors
- ◆ Sanitary flush diaphragm pressure transducer with digital readout
- ◆ Product inlet strainer
- ◆ Self-contained seal quench system
- ◆ IQ/OQ documentation & execution
- ◆ FAT, SAT, onsite start-up & operator training





Microfluidics International Corporation  
90 Glacier Drive, Suite 1000  
Westwood, MA 02090, USA

Tel: 617-969-5452 • 800-370-5452  
Fax: 617-965-1213

Email: [mixinginfo@idexcorp.com](mailto:mixinginfo@idexcorp.com)  
Web: [www.microfluidics-mpt.com](http://www.microfluidics-mpt.com)

---

**Industry-Leading Insights & News:**  
<https://www.microfluidics-mpt.com/blog>



Material Processing Technologies